

Dual chromatographic analysis of pesticide residues in eggplants: a green approach to food safety monitoring

Received: 25 May 2025

Accepted: 31 October 2025

Published online: 27 November 2025

Cite this article as: Hegazy A.M., Abdelrahman M.M., Kamel M.R. *et al.* Dual chromatographic analysis of pesticide residues in eggplants: a green approach to food safety monitoring. *Sci Rep* (2025). <https://doi.org/10.1038/s41598-025-27091-w>

Amira M. Hegazy, Maha M. Abdelrahman, Mariam R. Kamel & Adel M. Ahmed

We are providing an unedited version of this manuscript to give early access to its findings. Before final publication, the manuscript will undergo further editing. Please note there may be errors present which affect the content, and all legal disclaimers apply.

If this paper is publishing under a Transparent Peer Review model then Peer Review reports will publish with the final article.

Dual Chromatographic Analysis of Pesticide Residues in Eggplants: A Green Approach to Food Safety Monitoring

Amira M. Hegazy^{1*}, Maha M. Abdelrahman¹, Mariam R. Kamel², Adel M. Ahmed³

¹ Pharmaceutical Analytical Chemistry Department, Faculty of Pharmacy, Beni-Suef University, Alshaheed Shehata Ahmad Hegazy St., Beni-Suef, 62514, Egypt.

² Farshout Central Hospital, Qena, Egypt.

³ Department of Pharmaceutical Analytical Chemistry, Faculty of Pharmacy, South Valley University, Qena 83523, Egypt.

* corresponding author email: amira.hegazy@pharm.bsu.edu.eg

Abstract

The determination of pesticide residues is crucial for safeguarding human health, ensuring food safety, and contributing to environmental sustainability. Chlorfenapyr (CLF), abamectin (ABM), and fenpyroximate (FNP) are widely employed in eggplant (*Solanum melongena* L.) cultivation for pest control. This study presents the first reported simultaneous analysis of these pesticides using two validated chromatographic methods: high-performance thin-layer chromatography (HP-TLC, Method A) and reversed-phase high-performance liquid chromatography (RP-HPLC, Method B). Method A employed a mobile phase of methanol/chloroform/glacial acetic acid/triethylamine (7.00:2.50:0.50:0.10, v/v) with densitometric UV detection at 225 nm. The retention factors (*R_f*) for CLF, ABM, and FNP were 0.13, 0.34, and 0.85, respectively. Method B utilized a mobile phase of acetonitrile/5 mM ammonium acetate buffer (pH 4.0) (70:30, v/v) at a flow rate of 1.20 mL/min, with UV detection at 255 nm. The retention times for CLF, ABM, and FNP were 3.4, 6.1, and 8.4 minutes, respectively. Both methods were validated according to the key principles of ICH guideline Q2(R2) for the validation of analytical procedures, demonstrating robustness for assaying the ternary pesticide mixture in pure forms, commercial formulations (with excipients), and spiked eggplant samples. Both developed methods exhibited excellent detection and quantification capabilities. Method A demonstrated LODs ranging from 0.0006 to 0.0027 µg/band and LOQs from 0.0019 to 0.009 µg/band. Method B achieved even lower detection limits, with LODs from 0.0029 to 0.0296 µg/mL and corresponding LOQs from 0.0089 to 0.089 µg/mL for the target pesticides.

Key Words: HP-TLC; RP-HPLC; Chlorfenapyr; Abamectin; Fenpyroximate; eggplant.

INTRODUCTION

The Imperative for Pesticide Residue Monitoring

The environmental and health risks posed by pesticide residues have garnered significant global attention in both environmental sustainability and food safety research [1]. The potential for elevated pesticide levels remaining in food due to improper application has prompted governments and international organizations to implement strict maximum residue limits (MRLs) to protect consumer health.

Eggplant as a Key Crop and Its Vulnerability

Eggplant (*Solanum melongena* L.) is a regionally and commercially significant crop and one of the most popular traditional foods in Egypt [2], valued for its nutritional composition, including vitamins, proteins, and carbohydrates [3]. However, the susceptibility of eggplant fruits to insect pests often necessitates pesticide application by farmers to prevent substantial losses in productivity and quality [4]. Nevertheless, the extensive use of these pesticides can result in consumer exposure to potentially harmful residues [5, 6]. Therefore, monitoring residues prior to human consumption is critical, necessitating the development of efficient, cost-effective, and simultaneous analytical methods to determine residues of commonly used agrochemicals like chlorfenapyr (CLF), abamectin (ABM), and fenpyroximate (FNP).

Profile of the Target Pesticides

- 1- CLF is Chemically named [4-bromo-2-(4-chlorophenyl)-1-ethoxymethyl-5-(trifluoromethyl)pyrrole-3-carbonitrile], (Figure 1). It is a broad-spectrum pyrrole insecticide-acaricide that presents as a white-to-tan powder. It functions as a pro-insecticide, undergoing metabolic activation into a toxic compound upon pest ingestion, while demonstrating a relatively low impact on

- beneficial predatory insects [7]. It provides effective control against various pest species, including *Liriomyza* sp., *Frankliniella occidentalis*, *Spodoptera exigua*, and *Spodoptera litura* [8].
- 2- ABM is a macrocyclic lactone, (Figure 1), composed of a mixture of avermectin B1a (80%) and B1b (20%) homologs [9, 10]. It functions as a broad-spectrum, systemic insecticide that effectively controls mites and various mobile insect stages on fruits and vegetables [11]. Since its introduction in the 1980s, the use of ABM has significantly increased. Its primary variations, B1a and B1b, differ in their methylation patterns [12].
 - 3- FNP is an acaricide characterized by an oxime-bearing pyrazole structure [13, 14], (Figure 1). It exhibits high efficacy against larval stages by inhibiting mitochondrial electron transport [15]. Chemically designated as α -(4-phenoxyphenyl)- α , α -dimethyl-1H-pyrazole-3-propanenitrile, FNP is a pyrazole acaricide and insecticide classified within the sulfonanilide group [16, 17].

Analytical Gap and Literature Review

A review of existing literature reveals numerous chromatographic methods for the determination of CLF, ABM, and FNP individually or in combination with other pesticide classes. For instance, CLF has been determined via HPLC in various vegetables like cabbage [19, 20] and green beans [21], in formulations [22], and in water by GC [23]. ABM has been extensively estimated, including recent HPLC methods in 2024 for lake water, milk, and peach juice [24], and previously in avocados [25, 26], citrus fruits [27, 28], and other matrices [29-33]. Similarly, FNP has been determined by HPLC with other pesticides [34, 35] and in various samples including apples [36, 37], grapes [38], honey [39] and surface water [40].

Rationale and Novelty of the Present Study

Critically, no method has been reported for the simultaneous quantification of this specific ternary combination (CLF, ABM, and FNP). This represents a significant analytical gap, as these pesticides are commonly used in combination or rotation in eggplant cultivation to combat complex pest infestations and prevent resistance.

Consequently, current monitoring programs require laboratories to run multiple, separate single-analyte methods, which are time-consuming, resource-intensive, and costly. As pesticides are among the most toxic contaminants in fruits and vegetables [41], and with a growing trend towards multi-residue analysis [42-45], there is a clear need for a simultaneous method. This gap motivated our research team to develop two selective and sensitive chromatographic methods: high-performance thin-layer chromatography (HP-TLC, Method A) and reversed-phase high-performance liquid chromatography (RP-HPLC, Method B) for the estimation of this target mixture.

Strategic Application of Dual Chromatographic Methods

The strategic application of both techniques addresses different needs in the analytical workflow. HP-TLC serves as an efficient and economical preliminary screening tool with high throughput and minimal solvent consumption, making it ideal for initial monitoring. Conversely, RP-HPLC is preferred for its superior sensitivity, automation, and regulatory compliance for precise confirmatory quantification. This dual approach allows laboratories to maximize efficiency and resource utilization, using HP-TLC for initial screening followed by RP-HPLC confirmation for non-compliant samples, thereby balancing affordability with analytical rigor. This is particularly relevant in the current regulatory climate, where MRLs for such substances are under continuous review [46]. The newly developed methods were successfully applied to determine the investigated pesticides simultaneously in pure forms, commercial formulations, and spiked eggplant samples, fulfilling a critical need for monitoring this specific pesticide combination.

2. EXPERIMENTAL

2.1 Instruments

Method A

Data were acquired using a Camag TLC system that included a Linomat V auto-sampler and a TLC scanner III, operated with winCATS version 1.4.4.6337 software (CAMAG, Muttenz, Switzerland). The system

requirements included a deuterium lamp as the radiation source, absorbance mode for scanning, output of the chromatogram as an integrated peak area, slit dimensions of 4 x 0.45 mm, a Linomat V autosampler with a 100 μ L syringe for sample application, and a scanning speed of 20 mm/s.

Pre-coated 0.25-mm silica gel 60 F254 TLC aluminum plates (20 x 20 cm; Merck, Germany) were used as the stationary phase. A 4000-rpm electrical centrifuge (Zjmzym, China), a 0.1-100 μ L Rongtai variable-volume micropipette (Shanghai, China), and a 250 VM vortex mixer from Hwashin (Seoul, Korea) were used.

Method B

An Agilent Technologies 1260 Infinity series LC system connected to a 1260 UV-VIS detector was used, equipped with an Agilent Technologies Eclipse Plus C18 column (25 cm x 4.6 mm i.d., 5 μ m particle size).

Sample preparation

Other equipment included for sample preparation

- * Rotavapor R-300 (BUCHI Co., Switzerland),
- * Digital Ultrasonic Cleaner (USA),
- * Benchtop pH meter BP-3001 (TTBH Pte Ltd, Singapore),
- * Vortex mixer VSM-3 (Shelton Scientific, USA),
- * Thermo Scientific Megafuge-8 Benchtop Centrifuge (UK),
- * djustable micropipette with variable volumes (0.10–1000 μ L, Shandong, China), &
- * Vertical ultrafreezer ULF500 (-86 °C, Infrico medicare, Spain).

2.2 Material and Reagents

2.2.1 Pure samples:

The concerned pesticides (CP) were purchased from Sigma-Aldrich (Talat Harb, Cairo, Egypt). Their purities were certified to be 99.21%, 99.05%, and 99.10% for CLF, ABM, and FNP, respectively, according to the supplier's certificate of analysis. (Note: The text said "99.21, 99.05,

99.10 and 99.47 %" which lists four values for three pesticides. I corrected this to three values. Please verify the correct purities.)

2.2.2 Commercial formulations:

- Baylora® (Batch No. 502023), labeled to contain 24 % w/w CLF, was manufactured by Hebei Guanwang Agro Chemical, China.
- Aquachem® (Batch No. 2023115), labeled to contain 5.7 % w/w ABM, was manufactured by Jialingsu Subin Agro Chemical, China.
- Ortus Super® (Batch No. 2051947), labeled to contain 5 % w/w FNP, was manufactured by Shoura Chemicals, Cairo-Alexandria Desert Road, Egypt.

2.2.3 Chemicals and reagents:

Method A

Methanol (MeOH) and acetonitrile (ACN) were HPLC-grade reagents purchased from Sigma-Aldrich Chemie GmbH, Germany.

Method B

Acetone (ACT), chloroform (CFM), glacial acetic acid (GAA), ammonium acetate (AA), and triethylamine (TEA) were of analytical grade (El-Nasr Pharmaceutical Chemicals Co., Cairo, Egypt). De-ionized water (DW) was obtained from SEDICO Co., Egypt.

2.3 Standard Solutions and Sample Preparations:

2.3.1 Standard solutions:

Stock standard solutions of CLF, ABM, and FNP were prepared in MeOH at a concentration of 1000 µg/mL. Working standard solutions (100 µg/mL) were prepared by appropriate dilution. All stock standard solutions were freshly prepared on the day of analysis and stored in a refrigerator for no longer than 24 hours.

2.3.2 Field sampling and pesticide application

Field Trial Design: A controlled field trial was conducted in a representative eggplant (*Solanum melongena* L.) cultivation area in the El-Fayoum governorate, Egypt, during the winter season of 2024-2025. The trial plot area was 0.5 hectares.

Pesticide Application Regime: The commercial formulations Baylora® (CLF), Aquachem® (ABM), and Ortus Super® (FNP) were applied as a tank mixture at their recommended field rates (250 g/L for CLF, 50 mL/L for ABM, 100 mL/L for FNP) using a motorized knapsack sprayer. Application was performed at the fruiting stage according to standard agricultural practices. Environmental conditions during application were recorded: temperature was, $25 \pm 2^\circ\text{C}$, relative humidity was $60 \pm 5\%$, and wind speed was negligible.

Sample Collection: Sampling was performed according to the guidelines of the Joint FAO/WHO Meeting on Pesticide Residues (JMPR) [47]. Samples were collected at six intervals: 2 hours (0-day), and 1, 3, 7, 10, and 14 days after application (DAA). This allows for the construction of dissipation curves.

At each interval, three replicate samples were collected from randomly selected plants across the treated plot. Each replicate sample consisted of approximately 1 kg of fresh, ripe, marketable eggplants, yielding a total of ~3 kg per time point. To avoid cross-contamination, samples were handled with clean gloves. A control sample was collected from an untreated plot maintained under the same agronomic conditions before the application of pesticides.

Immediately after collection, samples were labeled, placed in sterile polyethylene bags, transported to the laboratory in an icebox (maintained at $\sim 4^\circ\text{C}$), and processed upon arrival.

2.3.3 Spiked eggplant (quality control) sample

To confirm the functional applicability of the calculated Limits of Quantification (LOQ) within the complex matrix, a single quality control (QC) sample was prepared by spiking the blank eggplant matrix extract with the three analytes at a concentration equal to the highest LOQ value determined by both methods.

2.3.4 Commercial formulations

An amount equivalent to 100 mg of CLF, ABM, and FNP from Baylora® powder, Aquachem®, and Ortus Super®, respectively, was accurately weighed and transferred separately to 100-mL volumetric flasks. Then, 75 mL of MeOH was added to each flask. The solutions were sonicated for 20 min, cooled, diluted to volume with MeOH to obtain 1 mg/mL stock solutions, and then filtered. Appropriate dilutions of the stock solutions were made to obtain 100 µg/mL working solutions.

2.4 Procedures

2.4.1 Chromatographic conditions

Method A

HP-TLC separation was performed on (20 x 10 cm) silica gel 60 F254 pre-coated plates (0.25 mm thickness). The plates were pre-washed with MeOH and activated at 100°C for 15 minutes. Samples were applied as 6-mm bands, spaced 5 mm apart and 10 mm from the plate's lower edge, using a CAMAG Linomat-V automated sampler.

Ascending linear development was performed in a chromatographic tank pre-saturated for 30 minutes with the mobile phase MeOH:CFM:ACT:GAA:TEA (7.00:2.50:0.50:0.10:0.10, v/v). Development proceeded to a distance of 9 cm at ambient temperature. After air drying, the plates were scanned using a CAMAG HP-TLC scanner at 255 nm. Data acquisition and analysis were performed using WINCATS software with a deuterium lamp. The peak areas were recorded.

Method B

Separation was performed using an Agilent Zorbax Eclipse Plus-C18 column (250 x 4.6 mm, 5µm). The mobile phase consisted of ACN:5 mM AA buffer (70:30, v/v), with the buffer pH adjusted to 4.00. The flow rate was 1.20 mL/min, and UV detection was performed at 255 nm.

Preparation of 5 mM AA buffer (1 Liter): 385.40 mg of AA (CH₃COONH₄, purity > 99%) was weighed into a beaker, dissolved in approximately 800 mL of DW, equilibrated at room temperature (20-

25°C), and adjusted to pH 4.0 with acetic acid before diluting to 1 L with DW.

2.4.2 Construction of calibration curves

Method A

Aliquots from standard stock solutions (1 mg/mL) of CLF (0.005–0.5 µg), ABM (0.002–0.1 µg), and FNP (0.01–0.12 µg) were accurately transferred to a set of 10-mL volumetric flasks and diluted to volume with MeOH. A 10 µL volume of each prepared solution was applied to HP-TLC plates (20 × 10 cm, pre-washed with MeOH and dried at 60°C for 5 min) as 6-mm bands using a Camag Linomat IV applicator. The bands were applied 10 mm from the bottom edge and 5 mm apart.

Linear ascending development was performed in a chamber pre-saturated for 30 minutes with MeOH:CFM:ACT:GAA:TEA (7.00:2.50:0.50:0.10:0.10, v/v) as the mobile phase. UV detection was performed at 255 nm.

Method B

Accurate aliquots of CLF, ABM, and FNP working standard solutions (100 µg/mL) were transferred into separate sets of volumetric flasks to prepare solutions with concentration ranges of 0.02–2.00, 0.009–1.00, and 0.09–1.20 µg/mL, respectively. Triplicate injections were performed for each concentration. The calibration curve was constructed by plotting the integrated peak area against the corresponding concentration, and regression equations were derived.

Chromatographic separation was performed in isocratic mode on a C18 column with a mobile phase of ACN: 5 mM AA buffer (70:30, v/v, pH 4.00). The flow rate was 1.20 mL/min, and UV detection was at 255 nm.

2.4.3 Repeatability

A minimum of nine determinations covering the specified procedure range (three concentrations/three replicates each) was investigated.

2.4.4 Reproducibility

Reproducibility was assessed through an inter-laboratory trial. Investigation of reproducibility is not typically required for regulatory submission but should be considered for the standardization of an analytical procedure, such as for inclusion in pharmacopoeias.

2.4.5 Field sample preparation:

The Eggplant peel extracts were prepared from samples collected in Egypt in Feb. 2025. Peels were manually removed, chopped, and a representative $15.0 \text{ g} \pm 0.1 \text{ g}$ portion of the homogenate was weighed into a 50-mL Teflon centrifuge tube then combined with 25 mL of CAN and centrifuged for 4 minutes. After adding 5 g of anhydrous sodium sulfate, the mixture was centrifuged at 6000 rpm for 14 minutes at 18°C. The resulting supernatant was collected and concentrated to 3 mL under vacuum at 40°C using a rotary evaporator. The concentrated extract was then filtered through a $0.45 \text{ }\mu\text{m}$ filter. Stringent laboratory practices, including the use of clean equipment, were employed throughout the preparation process to minimize cross-contamination. Samples were analyzed immediately after preparation.

2.4.6 Application to commercial formulations

The procedure described under the construction of the calibration curve was followed for the prepared solutions of the commercial formulations. The concentrations of the CP were calculated from the corresponding regression equations.

2.4.7 Application to field samples

The procedure described under the construction of the calibration curve was followed for the prepared eggplant samples. The concentrations of the CP were calculated from the corresponding regression equations.

To prevent cross-contamination, all equipment was meticulously cleaned between samples, and the analysis order included solvent blanks and control matrix samples to monitor for any carry-over, which was not observed.

3. Results and Discussion

3.1 Method development and optimization

Method A

Optimization of the chromatographic separation for the ternary pesticide mixture involved evaluating several mobile phase compositions. Initial trials with MeOH:CFM (1:9 and 4:6, v/v) and hexane:MeOH (8:2 and 4:6, v/v) yielded unsatisfactory resolution. A subsequent attempt using MeOH:CFM:acetic acid (7.00:2.50:0.20, v/v/v) resulted in moderate separation but exhibited peak tailing. Optimal separation was achieved with a mobile phase consisting of MeOH:CFM:ACT:GAA:TEA (7.00:2.50:0.50:0.10:0.10, v/v). Scanning at 255 nm was employed to maximize sensitivity. The chromatographic chamber was pre-saturated with the mobile phase for 30 minutes to ensure homogeneity and minimize volatility. Under these conditions, the R_f values for CLF, ABM, and FNP were determined to be 0.13, 0.34, and 0.85, respectively (Figure 2).

Method B

Several mobile phases were tested using an Agilent Zorbax Eclipse Plus-C18 (250 X 4.6 mm, 5 μ m) column. Starting with MeOH:DW (90:10, v/v), overlapping broad peaks for CLF and ABM were obtained. Adjusting the ratio of MeOH to DW from (80:20, v/v) to (70:30, v/v) and finally (60:40, v/v) did not resolve the overlapping CLF and ABM peaks. Using MeOH:0.02 M Na_2HPO_4 buffer (70:30, v/v, pH 4 with ortho-phosphoric acid) also failed to resolve the peaks and introduced significant noise in the chromatogram between 2.50 and 5.25 minutes. Replacing the phosphate buffer with 5 mM AA buffer improved the resolution but some noise remained. Replacing MeOH with ACN minimized the noise but further improvement was needed. Another system, ACN:5 mM AA buffer (90:10, v/v, pH 5.60), resulted in good separation but FNP eluted at 14.00 minutes and the ABM peak showed tailing. Finally, ACN:5 mM AA buffer (70:30, v/v, pH 4.00) provided good resolution and symmetric peaks for all analytes, with retention times of 3.4, 6.1, and 8.4 min for CLF, ABM, and FNP, respectively.

The flow rate was adjusted from 0.50 to 2.00 mL/min; complete separation within an acceptable time was achieved at 1.20 mL/min. Scanning was performed at various wavelengths (215, 220, 225, 230, 240, 255, 260, and 270 nm). A wavelength of 255 nm provided the best

sensitivity for the CP with the least noise, as depicted in Fig. 3. The run time was short, at less than 10 minutes.

3.2 Method validation

Both methods were validated according to the key principles of ICH guideline Q2(R2) [48] for the validation of analytical procedures. The validated parameters included linearity, accuracy, precision (repeatability and intermediate precision), specificity, LOD, LOQ, and robustness. It should be noted that reproducibility, as defined by an inter-laboratory trial, was not conducted as the primary focus of this work was on method development and initial validation.

3.2.1 Linearity

The linearity of an analytical method is crucial for ensuring the detector response is directly proportional to the analyte concentration across a defined range, allowing for accurate quantification. For reliable pesticide residue analysis, the linear range must cover all expected concentrations, from the LOQ to the highest anticipated residue levels. Crucially, the established MRLs for specific pesticides, such as those for CLF (0.821 mg/kg), ABM (0.05 mg/kg), and FNP (0.3 mg/kg), are officially declared and regularly updated by regulatory bodies like the Codex Alimentarius Commission and the European Union [49-51]. Therefore, it is imperative that the developed analytical method demonstrates excellent linearity across a range that encompasses these MRLs, as shown in Table 1.

3.2.2 Accuracy

Accuracy was assessed by calculating the percent recovery (%R) of known amounts of pure CP using the corresponding regression equations. Acceptable recovery percentages were obtained for both methods (Table 1). The standard addition technique was performed by adding known amounts of pure CP to their pharmaceutical formulations; good recoveries were obtained, showing no interference from excipients (Table 2).

3.2.3 Precision

Repeatability

Three different concentrations of pure pesticides were analyzed in triplicate (intra-day) using Methods A and B. Good RSD% values were obtained, verifying the methods' repeatability (Table 1).

Intermediate precision

The procedure was repeated inter-day (three consecutive days) for the analysis of the three chosen concentrations. Acceptable RSD% values were obtained (Table 1).

3.2.4 Specificity

Specificity of the methods was critically assessed by comparing chromatograms of:

- * Standard solutions of the individual pesticides and the ternary mixture.
- * Blank eggplant extract (control sample) to identify any potential co-eluting matrix interferences.
- * Eggplant extract spiked with the target pesticides at the LOQ and MRL levels.

As shown in Figure (2) (method A) and Figure (3) (method B), the blank matrix samples showed no significant interfering peaks at the retardation factors ($R_f=0.13, 0.34, \text{ and } 0.85$) or retention times ($R_t=3.4, 6.1, \text{ and } 8.4$) corresponding to CLF, ABM, or FNP.

Furthermore, the peaks for all three analytes in the spiked samples were well-resolved and symmetrical, confirming the methods' selectivity in the presence of the complex eggplant matrix. The high percentage recovery values (Table 1) further corroborate the lack of matrix-induced enhancement or suppression.

3.2.5 Limits of detection and quantitation (LOD and LOQ)

For Method A, LOD and LOQ were based on the standard deviation of the response and the slope of the calibration curves: $LOD = 3.3 \times SD/slope$; $LOQ = 10 \times SD/slope$. LOD and LOQ were also calculated

using a visual non-instrumental method according to ICH recommendations [48]. For Method B, the low LOD and LOQ values demonstrate the high sensitivity of the proposed method (Table 1).

3.2.6 Robustness

The robustness of the proposed methodologies was verified by introducing small, deliberate changes in the chromatographic conditions (e.g., change in acetic acid amount $\pm 0.02\%$, saturation time ± 5 min in Method A; change in organic modifier $\pm 1\%$, flow rate ± 0.05 mL/min in Method B). These changes did not significantly affect the R_f values, peak symmetry, or peak area (Table 1).

3.2.7 System Suitability

Overall system suitability testing was performed. Acceptable results for the system suitability parameters are shown in Table 1.

3.3 Application to commercial formulations

Methods A and B were applied to determine the CP in their commercial formulations. The acceptable recovery values verified the efficacy of the proposed methods for determining the pesticides in their commercial formulations (Table 2).

3.4 Field sample analysis

Field sample analysis showed initial concentrations of the studied pesticides (CLF, ABM, and FNP) to be 2.09, 0.06, and 0.08 mg Kg⁻¹ on the first day, which decreased to 0.082, 0.003, and 0.0108 mg Kg⁻¹ by day 10, respectively. These values are within the recommended doses for human health and vegetable productivity (Table S1).

Analysis of the blank control samples collected from the untreated plot confirmed the absence of any detectable residues of the target pesticides (concentrations below the LOD for all analytes), ensuring that the results from the treated plot were due solely to the applied pesticides and not background contamination.

3.5. Validation of LOQ in Eggplant Matrix

The calculated Limits of Quantification (LOQs) for all three analytes were functionally confirmed in the eggplant matrix to establish the lowest concentration that could be accurately and precisely quantified. A dedicated LOQ confirmation samples were prepared by spiking the blank eggplant matrix extract at a concentration equal to the highest calculated LOQ value among the three analytes. These spiked samples were analyzed in three replicates ($n=3$) alongside the other validation samples, by the developed methods. Table [3] presents the mean recoveries, standard deviations (SD), and relative standard deviations (RSD) for the quality control (QC) concentration with three replicate. This detailed data clearly validates the method's excellent reproducibility (all RSDs $< 2.5\%$) and accuracy (recovery 95–105%) within the complex eggplant matrix, directly supporting the feasibility of the claimed LOQs. This result, which easily meets the validation criteria for trace analysis, unequivocally demonstrates that the new methods are both accurate and reproducible at the quantification limit within the realistic matrix environment.

3.6 Statistical analysis

The results obtained from the proposed methods for the assay of pure samples of CLF, ABM, and FNP were statistically compared with those obtained from **reference HPLC methods** [7, 29, and 36].

A pair of hypothesis tests was employed for each compound to compare the accuracy and precision of the proposed methods against the reference method:

Student's t-test was used to compare the **means** (a test for bias/accuracy). A calculated t-value below the theoretical critical value (at $p = 0.05$) indicates that there is no statistically significant difference between the mean results obtained by the two methods.

F-test was used to compare the **variances** (a test for precision). A calculated F-value below the theoretical critical value (at $p = 0.05$)

indicates that there is no statistically significant difference in the precision of the two methods.

The calculated *t* and *F* values for all three analytes, as shown in Table 3, were found to be **below the theoretical critical values**. This provides statistical evidence that the proposed methods exhibit **equivalent accuracy and precision** to the established reference method, with no significant systematic error or difference in reproducibility. This statistical comparison is crucial as it validates the reliability of the new simultaneous method against a proven benchmark

3.7 Greenness considerations

3.7.1 A Comparative Evaluation of HP-TLC and RP-HPLC in Terms of Solvent Consumption and Waste

When considering the environmental impact of chromatographic methods for pesticide residue analysis, both (HP-TLC) and (RP-HPLC) present distinct advantages and disadvantages in terms of their "greenness." While RP-HPLC offers high resolution and sensitivity, its significant environmental drawback is its high consumption of organic solvents like acetonitrile and methanol, which generates large volumes of hazardous waste [52]. This necessitates efforts to "green" RP-HPLC by reducing column dimensions, flow rates, and run times, or by exploring less toxic "green solvents," though these can have specific limitations.

In contrast, HP-TLC generally offers a more environmentally friendly profile. It inherently uses significantly less mobile phase solvent than HPLC, as the solvent moves by capillary action on a thin layer plate, directly translating to reduced solvent consumption and waste generation. Furthermore, HP-TLC's ability to run multiple samples simultaneously on a single plate contributes to higher throughput and potentially lower energy consumption per sample. Therefore, for pesticide residue analysis, while RP-HPLC remains critical for its precision and scope, strategies like miniaturization and the adoption of HP-TLC-based approaches are increasingly valuable for mitigating the environmental footprint associated with routine pesticide residue analysis.

3.7.2 Sustainability Assessment of the Proposed Methods and Comparison with Reported HPLC Methods

To provide a quantitative assessment of these principles, the greenness of both proposed methods was evaluated using the Analytical GREEnness (AGREE) metric [53]. The AGREE calculator provides a comprehensive score based on the 12 principles of GREEN analytical chemistry, ranging from 0 (not green) to 1 (fully green).

The sustainability of the proposed analytical methods was comprehensively evaluated and benchmarked against existing HPLC methods by assessing their environmental impact and analytical performance. This evaluation was guided by green, red and blue chemistry principles, focusing on eco-friendliness, water conservation, and safety. Three modern assessment tools RAPI, BAGI, and AGREE were utilized to quantify and compare the methods' sustainability and analytical quality by Prajapati et al. [54 and 55].

3.7.2.1 Blueness Profile Evaluation

The proposed methods [7, 29, and 36] showed higher blueness scores (80 and 82.5) and more dark-blue segments in the BAGI pictograms, indicating better ecological alignment compared to the reported HPLC methods, which scored lower (57.5 to 67.5) with more light-blue segments (Table S2).

3.7.2.2 Greenness Profile Evaluation

The proposed methods obtained AGREE scores of 72 and 77, indicating superior environmental sustainability compared to the reported methods, which yielded lower scores of 58, 52, and 55 (Table S2).

3.7.2.3 Red Analytical Performance Index

Analytical quality and method robustness were additionally evaluated via the RAPI metric. The proposed methods achieved a high RAPI score of 80, signifying enhanced reliability and performance, surpassing the scores of previously reported methods (65, 57.5, and 60) (Table S2).

The proposed (Method A) achieved a score value which is higher than that of the proposed (Method B), confirming its strong environmental

profile due to minimal solvent use and low energy demands of it and reflecting the higher environmental impact associated with Method B's continuous solvent flow and higher energy consumption, despite an efficient mobile phase.

Therefore, for pesticide residue analysis, while RP-HPLC remains critical for its precision, the adoption of HP-TLC-based approaches is a valuable strategy for mitigating the environmental footprint. The choice of method, or their combination, effectively balances analytical performance with the principles of green analytical chemistry, as demonstrated by the AGREE assessment.

Conclusion

Selective, precise, and accurate HP-TLC and RP-HPLC methods were developed and validated for the simultaneous determination of chlorfenapyr (CLF), abamectin (ABM), and fenpyroximate (FNP). The HP-TLC method offers a key advantage: simultaneous analysis of multiple samples with minimal mobile phase consumption, enabling rapid and selective pesticide quantification. HP-TLC is cost-effective for pesticide screening (low solvent use, high throughput), while RP-HPLC offers higher sensitivity and regulatory compliance at a greater cost. Labs can combine both for optimal efficiency.

These chromatographic methods are suitable for local and cross-border inspection systems to monitor eggplant samples, ensuring compliance with pesticide residue limits for domestic use or export. The developed simultaneous methods provide a practical solution for regulatory laboratories tasked with monitoring this specific pesticide combination. Compared to running three separate single-analyte methods, our approach reduces analysis time, solvent consumption, labor, and overall cost per sample by approximately two-thirds. This enhances surveillance capabilities by enabling higher throughput screening of samples, ensuring more comprehensive compliance checking with MRLs for this common pesticide combination in a key crop like eggplant. Validation with commercial formulations confirmed that excipients do not interfere with pesticide quantification, demonstrating their utility for quality control of both bulk powders and formulated products.

Recommended Workflow for Regulatory Application

Accordingly, we outline a clear, step-by-step protocol:

1. Sample Preparation: Follow the extraction method detailed in Section 2.4.5.
2. High-Throughput Screening: Analyze all samples via the developed HP-TLC method (Method A).
3. Confirmation of Positive Hits: Any sample showing a peak near the MRL must be re-analyzed and confirmed using the more precise RP-HPLC method (Method B).
4. Reporting: Quantify using the HPLC method and report values.

Reference

1. Aktar MW, Sengupta D, Chowdhury A. Impact of pesticides use in agriculture: their benefits and hazards. *Interdiscip Toxicol*. 2009 Mar;2(1):1-12. doi: 10.2478/v10102-009-0001-7. PMID: 21217838;PMCID:PMC2984095.
<https://pmc.ncbi.nlm.nih.gov/articles/PMC2984095/>
2. Abdelbagi AO, Ismail REA, Ishag AESA, Hammad AMA. Pesticide Residues in Eggplant Fruit from Khartoum State, Sudan. *J Health Pollut*. 2020 Feb 28;10(25):200304. doi: 10.5696/2156-9614-10.25.200304. PMID: 32175175; PMCID: PMC7058141.
<https://pmc.ncbi.nlm.nih.gov/articles/PMC7058141/>
3. Rakha, M. Growth, yield, and fruit quality of Eggplant (*Solanum melongena* L.) as affected by irrigation intervals and foliar application of some antitranspirants. *J. Plant Prod*. 2014, 5, 2069–2083, DOI: 10.21608/jpp.2014.64846).
https://jpp.journals.ekb.eg/article_64846.html
4. Conn, K. Pepper and Eggplant Disease Guide Edited by Pepper & Eggplant Disease Guide; Seminis® Vegetable Seeds Inc.: Oxnard, CA, USA, 2006; pp. 1–72.
<https://doi.org/10.3390/agriculture13010116>
5. Malhat, F.; Abdallah, O.; Anagnostopoulos, C.; Hussien, M.; Purnama, I.; Helmy, R.M.A.; Soliman, H.; El-Hefny, D. Residue, dissipation, and dietary intake evaluation of FFNPyroximateacaricide in/on guava, orange, and eggplant under

- open field condition. *Front. Nutr.* 2022, 9, 1846.
<https://doi.org/10.3389/fnut.2022.939012><https://www.frontiersin.org/journals/nutrition/articles/10.3389/fnut.2022.939012/full>
6. Abdallah, O.I.; Malhat, F.M. Thiachloprid residues in green onion (*Allium cepa*) using micro liquid–liquid extraction and liquid chromatography–tandem mass spectrometry. *Agric. Res.* 2020, 9, 340–348. DOI: 10.1007/s40003-019-00440-8
 7. Cao Y, Chen J, Wang Y, Liang J, Chen L, Lu Y. HPLC–UV analysis of chlorfenapyr residues in cabbage and soil to study the dynamics of different formulations. *Sci Total Environ.* 2005 Nov 1;350(1-3):38-46. doi: 10.1016/j.scitotenv.2005.01.031. Epub 2005 Apr 21. PMID: 16227071.
<https://pubmed.ncbi.nlm.nih.gov/16227071/>
 8. Badawy MEI, Mahmoud MS, Khattab MM. Residues and dissipation kinetic of abamectin, chlorfenapyr and pyridaben acaricides in green beans (*Phaseolus vulgaris* L.) under field conditions using QuEChERS method and HPLC. *J Environ Sci Health B.* 2020;55(6):517-524.
doi:10.1080/03601234.2020.1726701. Epub 2020 Feb 17. PMID: 32065572. <https://pubmed.ncbi.nlm.nih.gov/32065572/>
 9. Campbell, W.C. Ivermectin and Abamectin; Springer Science & Business Media: Berlin/Heidelberg, Germany, 2012.
 10. The e-pesticide manual, British Crop Protection Council (2001) Wise and Loveys Information Services Ltd, Herts.
DOI:10.1365/s10337-007-0442-0.
<https://ouci.dntb.gov.ua/en/works/IDdkNNol/>
 11. Jansson, R.K., Dybas, R.A. (1998). Avermectins: Biochemical Mode of Action, Biological Activity and Agricultural Importance. In: Ishaaya, I., Degheele, D. (eds) *Insecticides with Novel Modes of Action*. Applied Agriculture. Springer, Berlin, Heidelberg.
https://doi.org/10.1007/978-3-662-03565-8_9
 12. R. P. Pohanish, *Sittig's Handbook of Pesticides and Agricultural Chemicals*, William Andrew, 2014.
 13. Dekeyser, M.A. Acaricide mode of action. *Pest Manag. Sci.* Former. *Pestic. Sci.* 2005, 61, 103–110
 14. Konno T, Kuriyama K, Hamaguchi H. “FFNPyroximate (NNI-850), a new acaricide,” in *Proceedings of the Brighton Crop Protection Conference, Pests and Diseases*. (Vol. 1), (Farnham: British Crop Protection Council) (1990). p. 71–8
 15. Naik R, Kale V, Dethe M. Studies on residues of flufenzin and fFNPyroximate on brinjal. *Int J Plant Protect.* (2009) 2:38–41,
<https://www.semanticscholar.org/paper/Studies-on-residues-of->

- [flufenzin-and-fFNPyrroximate-Naik-Kale/0d85d21097d8f2b2ebcfebbd25f701e63a50cdef](#)
16. National Center for Biotechnology Information (NCBI), PubChem Compound Summary for CID 9576412, FFNPyrroximate, 2025, retrieved February 18, 2025)
 17. European Food Safety Authority (EFSA), EFSA J., 2016, 14, 4382
 18. Xu F, Xu D, Hu M, Chen L, Xu C, Zha X. Dissipation behaviour, residue analysis, and dietary safety evaluation of chlorfenapyr on various vegetables in China. Food Addit Contam Part A Chem Anal Control Expo Risk Assess. 2022 Apr;39(4):724-739. doi: 10.1080/19440049.2021.2025269. Epub 2022 Feb 1. PMID: 35104200. <https://pubmed.ncbi.nlm.nih.gov/35104200/>
 19. Yongsong Cao, Jiuxin Chen, Yuelong Wang, Ji Liang, Lihua Chen, Yitong Lu HPLC/UV analysis of chlorfenapyr residues in cabbage and soil to study the dynamics of different formulations Science of The Total Environment Volume 350, Issues 1–3.2005 Pages 38-46, <https://doi.org/10.1016/j.scitotenv.2005.01.031>. <https://www.sciencedirect.com/science/article/abs/pii/S0048969705000720>
 20. Ahmed Gotah^{1,2}, Coralia V. Garcia¹, Sam-Pin Lee^{1,2}, Key Whang¹ An analytical method for quantifying diazinon, chlorfenapyr, and lufenuron in Napa cabbage using QuEChERS and HPLC-UV Food Science and Preservation- 2018;25(4):446-452- Jul 20, 2018. DOI: <https://doi.org/10.11002/kjfp.2018.25.4.446>https://www.ekosfop.or.kr/archive/view_article?pid=kjfp-25-4-446
 21. Badawy MEI, Mahmoud MS, Khattab MM. Residues and dissipation kinetic of abamectin, chlorfenapyr and pyridaben acaricides in green beans (*Phaseolus vulgaris* L.) under field conditions using QuEChERS method and HPLC. J Environ Sci Health B. 2020;55(6):517-524. doi: 10.1080/03601234.2020.1726701. Epub 2020 Feb 17. PMID: 32065572. <https://pubmed.ncbi.nlm.nih.gov/32065572/>
 22. Hafeez A, Tawab IA, Iqbal S. Development and Validation of an HPLC Method for the Simultaneous Determination of Fipronil, Chlorfenapyr, and Pyriproxyfen in Insecticide Formulations. J AOAC Int. 2016 Sep;99(5):1185-90. doi: 10.5740/jaoacint.16-0108. Epub 2016 Jul 26. PMID: 27464884. <https://pubmed.ncbi.nlm.nih.gov/27464884/>
 23. Wu, Y.; Zhou, Z. Application of liquid-phase microextraction and gas chromatography to the determination of chlorfenapyr in water samples. Microchim Acta 162, 161–165 (2008). DOI

- <https://doi.org/10.1007/s00604-007-0863-3>.
<https://link.springer.com/article/10.1007/s00604-007-0863-3#citeas>
24. Jie Zhao, Zi-Yu Li, Jia-Hui Zhu, Shi-Lei Ji, Yu-Shen Liu, Quan-Wen Liu, Lu-Liang Wang Determination of abamectin in food based on magnetic covalent organic frameworks and HPLC-FLD Journal of Food Composition and Analysis Volume 130- 2024- 106199. <https://doi.org/10.1016/j.jfca.2024.106199>.
<https://www.sciencedirect.com/science/article/abs/pii/S0889157524002333>
 25. Hernández Borges, J., Ravelo-Pérez, L.M., Hernández-Suárez, E.M. et al. Determination of Abamectin Residues in Avocados by Microwave-Assisted Extraction and HPLC with Fluorescence Detection. Chroma 67, 69–75 (2008). DOI: 10.1365/s10337-007-0442-0.
https://www.researchgate.net/publication/226472949_Determination_of_Abamectin_Residues_in_Avocados_by_Microwave-Assisted_Extraction_and_HPLC_with_Fluorescence_Detection
 26. Javier Hernández-Borges, Lidia M. Ravelo-Pérez, Estrella M. Hernández-Suárez, Aurelio Carnero, Miguel Ángel Rodríguez-Delga Analysis of abamectin residues in avocados by high-performance liquid chromatography with fluorescence detection Journal of Chromatography A, Volume 1165, Issues 1–2, 2007. <https://doi.org/10.1016/j.chroma.2007.07.033>.
<https://www.sciencedirect.com/science/article/abs/pii/S002196730701268X>
 27. Mohammad Rezaee, Hossein Ali Mashayekhi, Abolfazl Saleh, Yaser Abdollahzadeh, Mohammad Hosein Naeeni, Nazir Fattahi - Determination of abamectin in citrus fruits using SPE combined with dispersive liquid–liquid microextraction and HPLC–UV detection-Journal of SEPARATION SCIENCE- Volume36 , Issue 16 Pages: 2563-2708, August 2013, <https://doi.org/10.1002/jssc.201300356>.
<https://analyticalsciencejournals.onlinelibrary.wiley.com/doi/full/10.1002/jssc.201300356>
 28. A.I Valenzuela, M.J Redondo, Y Pico, G Font Determination of abamectin in citrus fruits by liquid chromatography–electrospray ionization mass spectrometry, Journal of Chromatography A, Volume 871, Issues 1–2, 2000 Pages 57-65. [https://doi.org/10.1016/S0021-9673\(99\)01190-5](https://doi.org/10.1016/S0021-9673(99)01190-5).
<https://www.sciencedirect.com/science/article/abs/pii/S0021967399011905>
 29. Xie, X., Wang, X. & Zhao, L. A Fast, Simple, and Reliable High-Performance Liquid Chromatography (HPLC) Method for

- Determining Abamectin Residues in Vegetables and Fruits. *Food Anal. Methods* 4, 203–211 (2011). DOI: 10.1007/s12161-010-9153-5.
https://www.researchgate.net/publication/225445528_A_Fast_Simple_and_Reliable_High-Performance_Liquid_Chromatography_HPLC_Method_for_Determining_Abamectin_Residues_in_Vegetables_and_Fruits
30. Xie X, Gong S, Wang X, Wu Y, Zhao L. Simplified RP-HPLC method for multi-residue analysis of abamectin, emamectin benzoate and ivermectin in rice. *Food Addit Contam Part A Chem Anal Control Expo Risk Assess.* 2011 Jan;28(1):19-25. doi: 10.1080/19440049.2010.527377. PMID: 21132591.
<https://pubmed.ncbi.nlm.nih.gov/21132591/>
 31. Yoshii K, Kaihara A, Tsumura Y, Ishimitsu S, Tonogai Y. Liquid chromatographic determination of emamectin, milbemectin, ivermectin and abamectin in crops and confirmation by liquid chromatography-mass spectrometry. *J Chromatogr A.* 2000 Oct 27;896(1-2):75-85. doi: 10.1016/S0021-9673(00)00466-0. PMID: 11093642. <https://pubmed.ncbi.nlm.nih.gov/11093642/>
 32. Sari, P., Razzak, M., & Tucker, I. G. (2004). Rapid, Simultaneous Determination of Levamisole and Abamectin in Liquid Formulations Using HPLC. *Journal of Liquid Chromatography & Related Technologies*, 27(2), 351–364. DOI: 10.1081/JLC-120027105.
https://www.researchgate.net/publication/233063751_Rapid_Simultaneous_Determination_of_Levamisole_and_Abamectin_in_Liquid_Formulations_Using_HPLC
 33. Bobby N. Brewer, Kevin L. Armbrust, Keith T. Mead, William E. Holmes Determination of abamectin in soil samples using high-performance liquid chromatography with tandem mass spectrometry-Rapid Communications in Mass Spectrometry Volume18, Issue15- 29 June 2004.
<https://doi.org/10.1002/rcm.1537>
<https://analyticalsciencejournals.onlinelibrary.wiley.com/doi/full/10.1002/rcm.1537>
 34. Faraneh Zareiyan, Habibollah Khajehsharifi, Hamid Reza Rajabi-Dispersive liquid–liquid microextraction based on deep eutectic solvent followed by high-performance liquid chromatography *Microchemical Journal* Volume 201,2024.
<https://doi.org/10.1016/j.microc.2024.110519>
<https://www.sciencedirect.com/science/article/abs/pii/S0026265X24006313>

35. Jongsung Noh, Seung-Woon Myung Determination of FFNPyroximate in Aqueous Matrices by Dispersive Liquid–Liquid Microextraction Bulletin of the Korean Chemical Society: Volume 41, Issue 1- January 2020. <https://doi.org/10.1002/bkcs.11922>.
<https://onlinelibrary.wiley.com/doi/full/10.1002/bkcs.11922>
36. Abd Al-Rahman, S.H., Almaz, M.M. & Osama, I.A. Determination of Degradation Rate of Acaricide FFNPyroximate in Apple, Citrus, and Grape by HPLC-DAD. Food Anal. Methods 5, 306–311 (2012). <https://doi.org/10.1007/S12161-011-9243-Z>.
https://www.academia.edu/109115672/Determination_of_Degradation_Rate_of_Acaricide_FFNPyroximate_in_Apple_Citrus_and_Grape_by_HPLC_DAD
37. B.L Halvorsen, C Thomsen, T Greibrokk, E Lundanes, Determination of fFNPyroximate in apples by supercritical fluid extraction and packed capillary liquid chromatography with Journal of Chromatography A Volume 880, Issues 1–2 2000. [https://doi.org/10.1016/S0021-9673\(00\)00293-4](https://doi.org/10.1016/S0021-9673(00)00293-4).
<https://www.sciencedirect.com/science/article/abs/pii/S002196730002934>
38. Malhat, F., El-Mesallamy, A., Assy, M., Madian, W., Loutfy, N. M., & Ahmed, M. T. (2013). Residues, half-life times, dissipation, and safety evaluation of the acaricide fFNPyroximate applied on grapes. Toxicological & Environmental Chemistry, 95(8), 1309–1317. DOI: 10.1080/02772248.2013.877245.
https://www.researchgate.net/publication/259601715_Residues_half_life_times_dissipation_and_safety_evaluation_of_the_acaricide_fFNPyroximate_applied_on_grapes
39. Jin Mun Kim, Seung-Woon Myung Determination of FFNPyroximate from Four Types of Honey by Liquid Chromatography/Tandem Mass Spectrometry Bulletin of the Korean Chemical Society: Volume 39, Issue 1 January 2018. <https://doi.org/10.1002/bkcs.11346>.
<https://onlinelibrary.wiley.com/doi/full/10.1002/bkcs.11346>
40. Geramizadegan, A., Ghazanfari, D., & Amiri, A. (2022). Determination of amount herbicide toxic fFNPyroximate in surface water by analysis molecularly imprinted solid phase extraction method and relative error assessment using artificial neural network model. International Journal of Environmental Analytical Chemistry, 103(20), 9002–9018. DOI: 10.1080/03067319.2021.2001465.
https://www.researchgate.net/publication/359286171_Determination_of_amount_herbicide_toxic_fFNPyroximate_in_surface_water_by_analysis_molecularly_imprinted_solid_phase_extraction_meth

[od and relative error assessment using artificial neural network mod](#)

41. Fenik, J., Tankiewicz, M., & Biziuk, M. (2011, June). Properties and determination of pesticides in fruits and vegetables. *TrAC - Trends in Analytical Chemistry*.
<https://doi.org/10.1016/j.trac.2011.02.008>.
42. Hegazy, A. M., Mahmoud, H. M., Elsayed, M. A., Ali, N. W., & Abdelfatah, R. M. (2023). Applying novel economic simple green sample preparation procedures on natural and industrial specimens for chromatographic determination of insecticidal residues. *Scientific Reports*, 13, 7209. <https://doi.org/10.1038/s41598-023-33421-7>
43. Zhao, P. et al. Rapid multipug filtration cleanup with multiple-walled carbon nanotubes and gas chromatography–triple–quadruple mass spectrometry detection for 186 pesticide residues in tomato and tomato products. *J. Agric. Food Chem.* 62, 3710–3725 (2014).
44. Wahab, S.; Muzammil, K.; Nasir, N.; Khan, M.S.; Ahmad, M.F.; Khalid, M.; Ahmad, W.; Dawria, A.; Reddy, L.K.V.; Busayli, A.M. Advancement and New Trends in Analysis of Pesticide Residues in Food: A Comprehensive Review. *Plants* 2022, 11, 1106. <https://doi.org/10.3390/plants11091106>.
45. Bairy, Partha. (2015). A Comparison Study of HPLC and HPTLC: Principles, Instrumentations and Applications. *ASIO Journal of Analytical Chemistry (ASIO-JAC)*. 1. 20-28.
46. Hegazy, A. M., Abdelfatah, R. M., Mahmoud, H. M., & Elsayed, M. A. (2020). Development and validation of two robust simple chromatographic methods for estimation of tomatoes specific pesticides' residues for safety monitoring prior to food processing line and evaluation of local samples. *Food Chemistry*, 306. <https://doi.org/10.1016/j.foodchem.2019.125640>.
47. FAO/WHO, "Recommended Methods of Sampling for Determination of Pesticide Residues," (1986) 8: 2nd edition
48. ICH: Harmonized Tripartite Guideline, "Validation of Analytical Procedures: Text and Methodology" Q2 (R1), 2005.
49. Codex Alimentarius Commission. Pesticide residues in food and feed. Food and Agriculture Organization of the United Nations. Accessed June 18, 2025. <https://www.fao.org/fao-who-codexalimentarius/codex-texts/dbs/pestres/en/>.
50. European Commission. EU pesticides database. Accessed June 18, 2025. https://food.ec.europa.eu/plants/pesticides/eu-pesticides-database_en.

51. Horticulture Innovation Australia. *Eggplant Strategic Agri-Investment Research Plan (SARP) 2021*. Hort Innovation; 2021.
52. Perrucci M, De Laurenzi V, Dainese E, Locatelli M, Ulusoy HI, Kabir A, et al. Green chemistry and green analytical chemistry. In: Locatelli M, ed. *Green Analytical Chemistry*. Elsevier; 2025:1-14.
53. Locatelli M, Kabir A, Perrucci M, Ulusoy S, Ulusoy HI, Ali I. Green profile tools: Current status and future perspectives. *Adv Sample Prep*. 2023;6:100068. doi:10.1016/j.sampre.2023.100068.
54. Pintu Prajapati, Bageshree Rana, Veera Shakar Pulusu, Shailesh Shah, Simultaneous Chromatographic Estimation of Vildagliptin and Dapagliflozin Using Hybrid Principles of White Analytical Chemistry and Analytical Quality by Design, *Journal of AOAC INTERNATIONAL*, Volume 107, Issue 1, January-February 2024, Pages 212–222, <https://doi.org/10.1093/jaoacint/qsad108>.
55. Prajapati, P., Rana, B., Pulusu, V. S., & Shah, S. (2023). Method operable design region for robust RP-HPLC analysis of pioglitazone hydrochloride and teneligliptin hydrobromide hydrate: Incorporating hybrid principles of white analytical chemistry and design of experiments. *Future Journal of Pharmaceutical Sciences*, 9(1), 1-15. <https://doi.org/10.1186/s43094-023-00546-5>.

Acknowledgment

The authors would like to thank for Dr. Mohamed El-Sayed for his effort in this work progress.

Credit author statements

A.H.: Project administration, editing, and reviewing the abstract and introduction parts as well as the article, investigation, and supervision. **M.A.:** Editing and reviewing the whole experimental part. **M.K.:** Resources, investigation, project funding, and writing (original draft). **A.A.:** Editing and reviewing the figures, tables, and supervision.

Funding

The Egyptian Ministry of Higher Education in cooperation with The Egyptian Knowledge Bank (EKB) was the base of the work funding.

Competing interests

The authors declare no competing interests.

Data availability

All data generated or analysed during this study are included in this published article.

ARTICLE IN PRESS

Figures

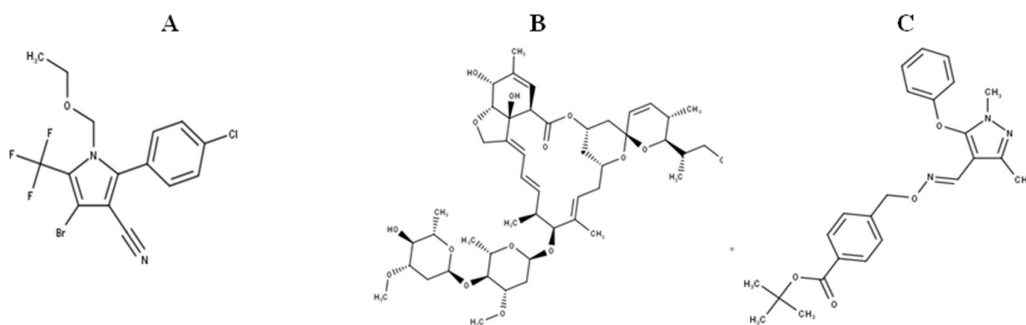


Figure 1. Chemical structures of **A**: CLF, **B**: ABM, and **C**: FNP.

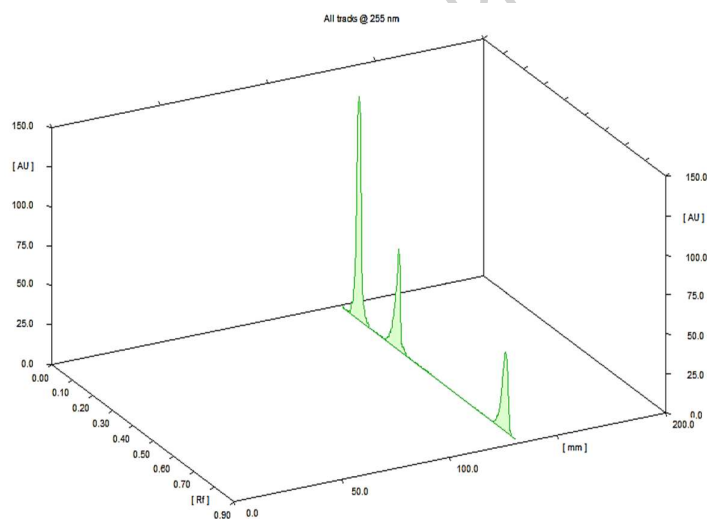


Figure 2. 3D HP-TLC-densitogram of the mixture: CLF (R_f=0.13), ABM (R_f=0.34), and FNP (R_f=0.85) using MOL: CM: AN: GAA: TEN (7.00: 2.50: 0.50: 0.1, v/v) as a MP and UV detection at 255 nm.

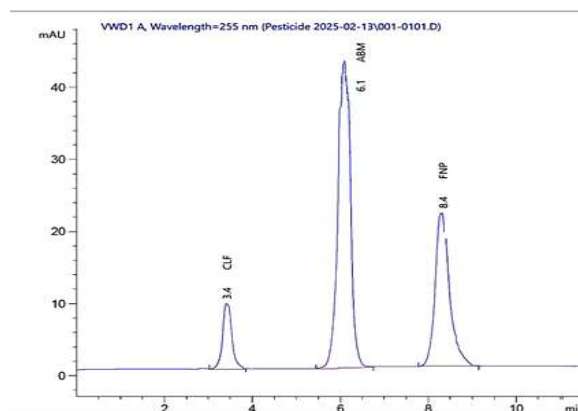


Figure 3: RP-HPLC chromatogram of the mixture: 0.40 $\mu\text{g/mL}$ CLF, 1.00 $\mu\text{g/mL}$ ABM, and 1.20 $\mu\text{g/mL}$ FNP using AN: 5 mM AA buffer (70: 30, v/v), pH adjusted to 4.00, Flow rate 1.20 mL/min, and UV detection at 255 nm.

Tables

Table1: Validation results and system suitability testing parameters of the methods A and B.

Validation criteria	Values measured						Approved criteria that follow guidelines of the ICH
	Method A			Method B			
	CLF	ABM	FNP	CLF	ABM	FNP	
Linearity,R=	0.99995	0.99997	0.99991	0.99997	0.99993	0.99992	r ≥ 099
Slope	2282.90	7439.49	5839.80	420.98	868.59	626.93	
Intercept	0.17	2.59	2.41	-1.79	1.33	-0.68	
Range	0.005-0.50	0.002-0.10	0.01- 0.12	0.02-2.00	0.0091.00	0.09-1.20	---
	(µg/band)			(µg/mL)			
Precision (6 replicates)	0.34	0.62	0.59	0.21	0.84	0.71	RSD ≤ 2%
Accuracy (Recovery means±SD)	100.24% ±1.043	100.52% ±1.014	100.19% ±0.754	100.17% ±0.702	99.77% ±1.159	99.99% ±0.998	100±2%
Specificity R _s	1.62	5.08		2.08	2.08		complete resolved peaks Resolution is > 1.5 not significantly interfere with quantification Using the formula: 3.3 x SD / slope Using the formula: 10 x SD / slope Every alteration should have a pooled RSD of less than 3%
ME %		-8.5			+6.2		
LOD	0.0015	0.0006	0.0027	0.0063	0.0029	0.0296	
		(µg/band)			(µg/mL)		
LOQ	0.0049	0.0019	0.009	0.0191	0.0089	0.089	
		(µg/band)			(µg/mL)		
Robustness	0.79	0.64	0.91	0.48	0.97	0.88	
System Suitability Testing Parameters	Method A			Method B			Reference value
	FNP	ABM	CLF	FNP	ABM	CLF	
Tailing factor* (T)	1.31	1.22	1.09	1.17	1.18	1.08	< 1.5
Capacity factor (K')	7.41	5.14	2.48	6.98	3.87	2.13	1- 10
Resolution (R _s)	2.71	3.01	5.08	1.62			> 1.5
Selectivity (α)	1.44	2.08	1.80	1.82			> 1
Column efficiency (N)	2254	1254	870	-	-		↑with efficiency of the separation
HETP (cm plate ⁻¹)	0.021	0.015	0.087	-	-		↓ the value↑ the column efficiency

HETP = height equivalent to theoretical plate, (cm.plate^{-1}).

*Calculated using three peaks.

Table 2: Determination of CLF, ABM and FNP in their commercial formulations Baylora®, Aquachem® Powder, Ortus super®, respectively, by the proposed methods and application of the standard addition technique.

Method A					Method B				
Taken (µg/band)	Found* %	Pure added	Pure Found	R %	Taken (µg mL ⁻¹)	Found* %	Pure added	Pure Found	R %
		(µg/band)						(µg/band)	
CLF									
0.04	99.87	0.03	0.03	100.01 99.98 100.12	0.40	100.08	0.30	0.29	99.98 98.01 100.54
		0.05	0.04				0.40	0.38	
		0.07	0.07				0.50	0.50	
		Mean± SD					100.05 ± 0.99	Mean± SD	
ABM									
0.02	99.87	0.03	0.02	99.78 100.21 99.81	0.10	99.69	0.30	0.30	100.12 100.01 98.24
		0.05	0.05				0.40	0.40	
		0.07	0.07				0.50	0.48	
		Mean± SD					99.54 ± 0.71	Mean± SD	
FNP									
0.06	100.34	0.03	0.03	99.62 99.81 98.69	0.40	100.25	0.30	0.31	100.99 99.67 100.89
		0.05	0.05				0.40	0.39	
		0.07	0.06				0.50	0.51	
		Mean± SD					99.10 ± 0.51	Mean± SD	

*Average of three determinations.

Table 3: Statistical comparison of the results obtained by the proposed methods and the reported methods for determination of pure CLF, ABM, and FNP.

Parameter	Method A			Method B			Reported Methods (7, 29 & 36)		
	CLF	ABM	FNP	CLF	ABM	FNP	CLF	ABM	FNP
Mean	100.24	100.52	100.19	100.17	99.77	99.99	99.38	99.20	99.41
SD	1.04	1.01	0.75	0.70	1.16	0.99	1.31	1.18	1.30
Variance	1.08	1.02	0.56	0.49	0.45	0.98	1.72	1.39	1.69
N	7	7	7	7	7	7	7	7	7
Student's <i>t</i>-test* (2.45)	0.57	0.87	0.45	0.12	0.89	0.09	----	----	----
<i>F</i>-test* (4.28)	1.98	1.68	2.14	1.89	2.01	1.68	----	----	----

*Figures in parenthesis are the corresponding tabulated values at $p = 0.05$.